

QUICKIE®  
IRIS

INTELLIGENT ROTATION IN SPACE



## EASY MANEUVERABILITY

(IRIS) Intelligent Rotation in Space technology (IRIS) means the seat frame rotates around the rider's center of gravity, resulting in a smooth tilting motion that is easy to control and allows for the shortest possible wheelbase for easy maneuverability.

## UP TO 55° OF TILT IN SPACE

A 55-degree tilt facilitates gravitational changes, enhancing sitting tolerance, digestion, and pressure distribution. The IRIS enables safe transfers, aids in repositioning hips and trunks, manages edema and blood pressure, and optimizes skin integrity by redistributing pressure.

## BUILT-IN ADJUSTABILITY AND GROWTH

Independent front & rear center of gravity and center of rotation adjustments can be made to optimize tilt performance and maneuverability when personal or additional medical items like an aftermarket backrest, custom mold, or ventilator are added.

## BROADEST RANGE OF OPTIONS TO SUIT COMPLEX NEEDS


Contracture footplate, dynamic back, recline, power tilt, and transit securements are just a few of the options available. The IRIS also offers the broadest range of lower extremity features, from enhanced hanger latches with increased durability for high-tone riders to swing-in/out hangers.

## PROVIDES THE POSTURAL SUPPORT REQUIRED

From uncomplicated to complex postural needs, bundle with JAY seating and positioning supports, RIDE Designs custom seating options, and a Whitmyer head support to customize optimal support even while rotating.

## POWER TILT

With the optional Power Tilt, the rider can effortlessly adjust their body position for maximum comfort, relaxation, and skin protection. With the Attendant Control, caregivers can tilt the rider without any manual effort.

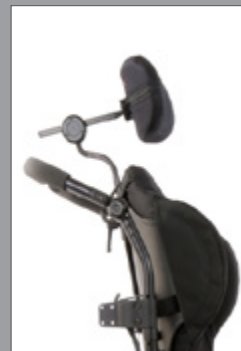

## TECHNICAL DATA

HCPCS CODE	E1161 / E1233
SEAT WIDTH	14" to 22"
SEAT DEPTH	14" to 22"
TILT RANGE	40° / 55° / 60°
SEAT HEIGHT FRONT	12.5" to 19.5"
SEAT HEIGHT REAR	12.5" to 19.5"
BACK HEIGHT ADJUSTMENT	15" to 21" or 18" to 24"
BACK ANGLE ADJUSTMENT	-30° to 5°
CHAIR WEIGHT	38 lbs. (w/foot supports, 18 x 18)
WEIGHT CAPACITY	300 lbs. / HD 400 lbs.
TRANSIT	Yes *

Refer to the order form as specifications may vary based on configuration.

\* Conforms with ANSI/RESNA WC-19 standards. Please refer to QUICKIE & ZIPPIE WC-19 Transit Securement Supplement (MK-100072) for additional information. A supplement is available for download at [www.SunriseMedical.com](http://www.SunriseMedical.com).

Learn more



 This document can be viewed in PDF format at [www.SunriseMedical.com](http://www.SunriseMedical.com).