

EVEN STRONGER STILL THE LIGHTEST

The Quickie Q7 NextGEN exemplifies our pursuit to optimize strength and use the lightest weight materials so you can experience ultimate performance day in and day out.





TOP OF THE LINE MATERIALS AND TECHNOLOGY

The genius behind the Q7's featherweight frame is 7000 series aerospace aluminum. ShapeLoc™ Technology— the process of heat treating 7000 Series aluminum— increases its strength, allows for thinner tube walls, and decreases overall weight. By combining ShapeLoc Technology with ovalized tubing, the Q7 delivers an incredibly strong, rigid ride that **weighs as little as 13.3 pounds with wheels** and maximizes your speed and efficiency with every push.



BALANCED PERFECTION

The Q7 NextGEN's performance geometry is designed for precision and grace. With custom frame lengths and caster wheel placement options, the Q7 properly distributes your weight over the rear wheels so you can be quick and nimble while maintaining stability.





Ride in Style

The Q7 NextGEN's minimal, ultra-modern design takes backstage and ensures YOU are in the limelight. Pick from 6 anodized caster color options and 30 frame colors including the Q7 NextGEN exclusive Glossy White.



Active Rigid Performance Frame Option

Ideal for those who know exactly what they need. A fully rigid backrest and axle plate reduce the overall weight by an average of one pound and transfer energy directly from your arms to the axles for the most efficient propulsion.

FINE TUNE YOUR
Q7 NEXTGEN
TO YOUR
INDIVIDUAL
ERGONOMIC AND
PERFORMANCE
REQUIREMENTS
AND YOU'LL BE
LOVING
LIFE IN THE
FAST LANE



Angle-Adjustable, Fold-down Backrest

Strong and easy to adjust from -27° to 8° to accommodate changing postural needs.




Caster Housing Design

Ovalized tubing, stylish hydroformed caster arm and a single-sided caster fork optimize balance, make the chair easier to turn and look really cool!



Xtender™ Power Assist Option

Developed by Quickie and Yamaha, the Xtender features motorized rear wheels that increase the force you apply to the handrims by four times!

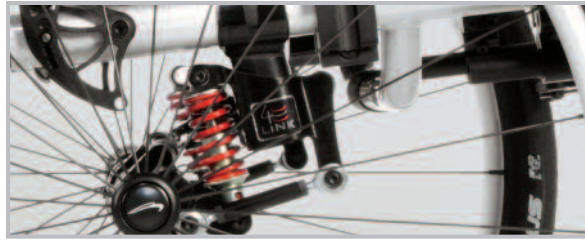

4-Link Rear Suspension

Sporty and lightweight suspension for a comfortable ride.

Q7 NEXTGEN SPECS

HCPCS Code: K0005

Lowest Overall Weight

Adjustable: 14.1 lbs
Active Rigid: 13.3 lbs

Weight Capacity

265 lbs

Seat Width

12" to 20"

Seat Depth

12" to 20"

Center of Gravity Adjustment

0" to 4"

Back Height Adjustment

8" to 20"

Front Seat Height

16" to 21"

Rear Seat Height

13" to 20"

Frame Angle

75°, 80°, 85°

Warranty

Limited/Lifetime