



J3[®] CUSTOM BACK SUPPORT

For those individuals who require a more complex seating solution, the JAY[®] J3 Custom Back Support can be used when an off-the-shelf back cannot accommodate an individual's orthopedic asymmetries. Featuring the J3 Back that you know and love, this mix-and-pour custom-contouring system makes it easy to create a personalized solution to suit your client's specific needs.



SUPERIOR CLINICAL SEATING

J3[®] CUSTOM BACK SUPPORT

Benefits of a J3 Custom Back Support

- Provides proper body positioning customized to the client's unique skeletal structure
- Utilizes JAY J3 shell, cover, and hardware
 - 3DX Microclimatic outer cover dissipates heat and moisture, increasing comfort and sitting tolerance
 - The zipper is conveniently located at the top of the back support for easier access during the mold-in-place process
- Coded as an E2617 for the custom seating market
 - Covered when a prefabricated seating system is not sufficient to meet seating and positioning needs

Fast and easy to create

- Mix the foam-in-place and pour into the molding bag within the J3 back custom cover secured on the client's chair
- The client leans into the bag and the foam will mold to the client's optimal seating position
- Available in three foam densities: extra soft, soft and medium

JAY TECHNOLOGY IN THE J3 CUSTOM BACKREST




POSITIONING




CONTOURING



FOUR WAY
STRETCH

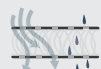

MICROCLIMATIC
MATERIAL

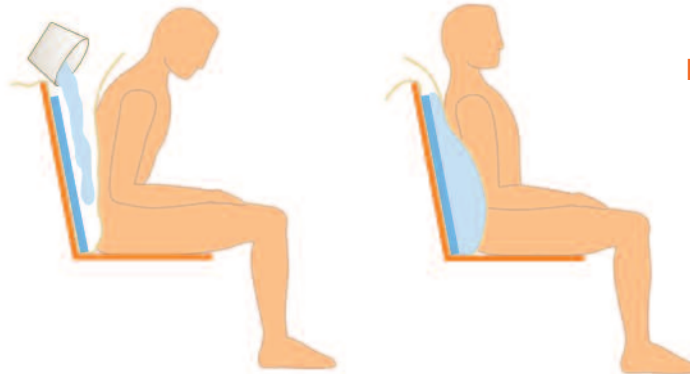

Foam-In-Place

The J3 Custom Back Support utilizes foam-in-place to create a personalized mold for optimal contact with the individual



Specifications:

HCPCS Code
E2617

Weight Capacity
300 lbs.

Width
Fits 14" to 20"

Contour
2.2"

Height
12" to 22.5"

Cover Material
3DX™ Spacer Fabric
with Lycra® Stretch

TH
Contact: 2-Point
Angle: +/- 20°
Width: 2"
Depth: 2.25"

FS
Contact: 4-Point
Angle: +/- 10°
Width: 2"
Depth: 1.38"

JAY Your Way
Customizations Available



Sunrise Medical (US) LLC, Fresno, CA 93727
MK-129738 Rev. A ©11.2019

800.333.4000 www.SunriseMedical.com

