



QUICKIE[®]

EXPERIENCE LIFE
FROM NEW HEIGHTS

Q700-UPTM

QUICKIE Q700[®]-UP M

with SEDEO[®] ERGO[™] SEATING SYSTEM

For those who want to experience life from new heights, discover Q700-UP M.

QUICKIE combined all the benefits of the SEDEO ERGO seating system with the power to stand whenever you want.



Q700 M MID-WHEEL DRIVE PERFORMANCE

The Q700 M mid-wheel drive base sets the bar for powered wheelchair technology, inside and out. Travel further with its superior drive performance. Feel secure over challenging terrain with SpiderTrac[®] 2.0 suspension and its independent movement of all 6 wheels. It's a driving experience like never before.

BENEFITS OF STANDING*

The Q700-UP M offers many health and lifestyle benefits:

- Improves pressure relief for those who may have limited ability to independently weight shift or who sit for long periods
- Enhances function and ability to access your environment
- Improves range of motion and reduces the risk of contractures
- Supports vital organ capacity and function and improves circulation
- Promotes bone health and reduces the occurrence of skeletal deformities
- Reduces abnormal muscle tone and spasticity

90° KNEE ANGLE


The centermount footrest gives a true 90° knee angle in seated position. The ergonomic knee pad is also depth adjustable.

* See NCART Standing Device Funding Guide available on www.sunrisemedical.com/EIM

ADJUSTABLE TO FIT YOU

The premium headrest (1) can be easily adjusted or removed. All seat widths and depths (2) are infinitely adjustable with just one tool. Also, the seat and backrest width (3) can be adjusted independently of each other. Even the armrest height and angle (4) are simple to adjust.

- Pressure Relief Reminders
- Programmable Memory Seating
- Footplate Options

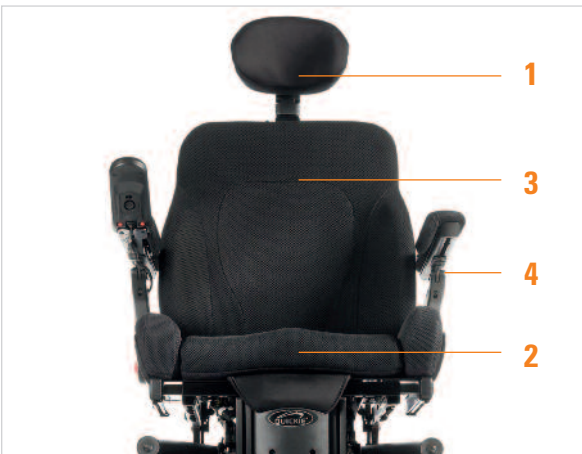

BIOMETRIC REPOSITIONING

Moving out of a reclined position often requires use of the tilt function to regain your proper seated posture. With the patented Biometric Repositioning, when you move from reclined to seated, the SEDEO ERGO will first pass through an intermediary tilted position, using gravity to realign your pelvis in one smooth, effortless motion.




ANTI-SHEAR BACKREST

Transitioning in and out of recline doesn't need to be uncomfortable. The SEDEO ERGO's power recline is supported by 7" of anti-shear to better maintain the positioning of the headrest, and even chin control, during recline.

MEMORY POSITIONS

With memory seating, you can save up to six programmable positions and transition to them with a simple press of a button.

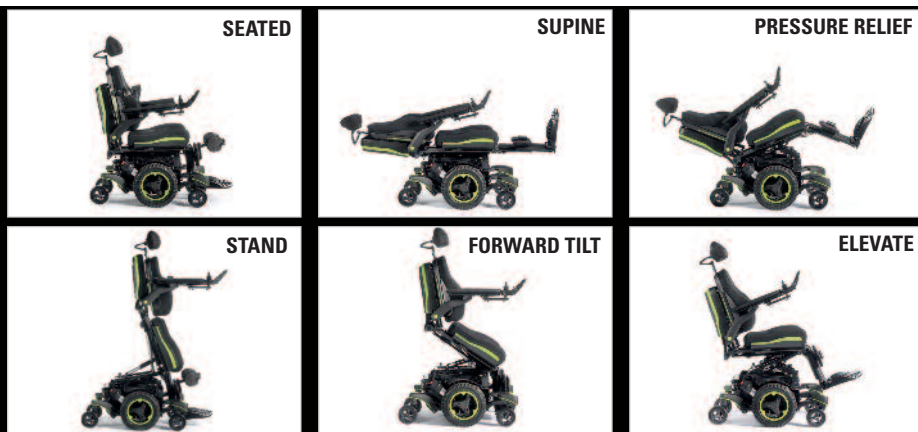

DISTINGUISHED STYLE

Choose from eight different shroud colors to reflect your personal style.

The Matte Black shroud also comes with a choice of six contemporary color inlays.



Q700-UP M	MPC	MPD
HCPC code	K0861	K0884
Max. user weight	265 lbs	
Seat width	16" to 20"	
Seat depth	17" to 22"	
Seat-to-floor height	18" to 19"	
Backrest height	22" to 26"	
Total width	24.5"	
Total length	36.5"	
Turning radius	19.8" <i>(excludes footplate)</i>	
Seat memory functions	Seated, Pressure Relief, Supine, Forward Tilt, Elevate, Stand	
Max. speed	6 mph <i>(standard)</i> , 8 mph <i>(optional)</i>	
Battery size	Group 24	
Overall weight	335 lbs <i>(w/o batteries)</i> , 435 lbs <i>(w/batteries)</i>	
Drive wheel position	Mid	
Ground clearance	3"	
Wheel options	14" drive wheels; 6" casters	
Max. range	up to 18 miles <i>(w/Group 24 batteries and 6 mph motors)</i>	
Transit	ISO 7176-19:2008 Transit Approved	
Warranty	Lifetime Limited	